



don whitley
scientific
excellence in microbiology

Whitley Workstation Single Plate Entry Systems



**Introduce individual samples
without compromising internal
atmospheric conditions**


**Versions for the majority of our
workstations can be retro fitted**

The Whitley Workstation range includes a variety of options and accessories enabling you to tailor a system to your needs and your laboratory working practices.

A very useful option is a single plate entry system, which allows individual Petri dishes to be introduced into the chamber without compromising internal atmospheric conditions. Multiwell plates, cell culture flasks and similar sized objects can be introduced using the same system.

The single plate entry system can be extremely handy for on-call and emergency use where individual samples can be introduced quickly without the need to use the airlock or ports.

More convenient working methods and lower gas usage, what more could you want?

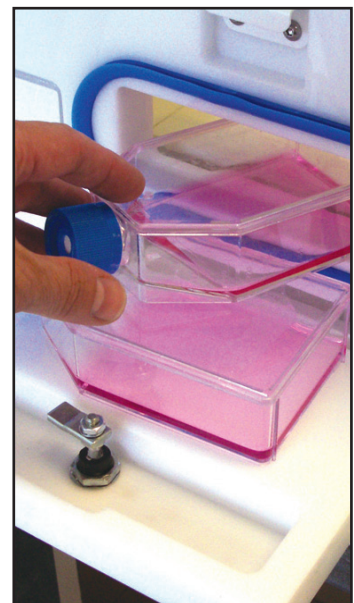

9cm - accommodates 90mm Petri dishes and items of a similar size						15cm	
	A02711	A03125	A03225	A03325	A07208	A06103	
DG250	-	●	●	●	-	-	DG250
A35	-	-	-	-	●	●	A35
A45	-	-	-	-	●	●	A45
A55	-	-	-	-	●	-	A55
A85	-	●	●	●	-	-	A85
A95	-	●	●	●	-	-	A95
VA500	●	-	-	-	-	-	VA500
VA1000	●	-	-	-	-	-	VA1000
H35	-	-	-	-	●	●	H35
H45	-	-	-	-	●	●	H45
H85	-	●	●	●	-	-	H85
	w: 98mm h: 20mm Front loading Factory fitted option	w: 107mm h: 26mm Front loading Can be retro fitted	w: 107mm h: 26mm Front loading With longer tray for larger items Can be retro fitted	w: 98mm h: 20mm Front loading With folding tray for an unobstructed work area Can be retro fitted	w: 107mm h: 26mm Front loading Can be retro fitted on units with a removable front	w: 152mm h: 52mm Located on the left hand panel of the unit Can be retro fitted	
Key:	- Not available			● Option available			

Single Plate Entry Systems are designed for occasional use and not as a direct alternative to airlocks or porthole systems that allow sample transfer.

In the interests of a policy of continuous product improvement the company reserves the right to alter product specifications without prior notice. All rights reserved. © 2014 Don Whitley Scientific Limited.

Don Whitley Scientific Limited

14 Otley Road, Shipley, West Yorkshire, BD17 7SE, England.

Tel: +44 (0)1274 595728 Fax: +44 (0)1274 531197

Website: www.dwsscientific.co.uk Email: info@dwsscientific.co.uk